

→ 限時特價，買一送一!!!

即日起至2025/11/30日止，可混搭(以低價者作為贈送)

Western blot 相關 Buffer

貨號	品名	包裝	價格
BR120-1L	10X TBS Buffer, pH7.6	1L	800
BR140-1L	10X Western Transfer Buffer	1L	800
BR150-1L	10X Tris Glycine SDS Running Buffer	1L	800
BR300-1L	4X Resolving Gel Buffer	1L	1,500
BR310-1L	4X Stacking Gel Buffer	1L	1,500
BR620-1L	10X MOPS SDS Running Buffer	1L	1,500
BR510-1L	10X TBST Buffer	1L	1,000
BR540-1L	10X PBST Buffer	1L	1,000
BS-20-500	Western Blot Stripping Buffer	500ml	2,700
BR590-500	Ponceau S Staining Solution	500ml	1,500
BR600-25	5X Protein Sample Buffer	25ml	2,500

註：含 Glycine 產品 (如TG buffer 或 TGS buffer) 若有變黃，屬於正常現象，請安心使用或置於陰暗處避光!

核酸實驗及相關 Buffer

貨號	品名	包裝	價格
BR410-100G/ 500G	Low EEO Agarose	100G/ 500G	4,000/ 9,000
BR160-1L	20X SSC Buffer, pH7.0	1L	800
BR170-1L	10X TE Buffer	1L	800
BR190-1L	10X TBE Buffer	1L	1,000
BR200-1L	50X TAE Buffer	1L	1,000
BR350-1L	20X SSPE Buffer 1L	1L	1,500
BR610-1L	10X MOPS Buffer	1L	1,500

IHC 相關產品

貨號	品名	包裝	價格
BR520-1L	10X Sodium Citrate Buffer, pH6.0 (Antigen retrieval buffer)	1L	1500
BR530-1L	10X Tris-EDTA Buffer, pH9.0 (Antigen retrieval buffer)	1L	1500
BR540-1L	10X PBST Buffer	1L	1,000
PHC02-20	I.S Circle Writer (liquid blocker) IHC 阻水筆	1支	2,000

Protein marker/ DNA Ladder 電泳標準品

貨號	品名	包裝	價格
BR331-500	100bp plus DNA RTU Ladder, 100 bp ~1.5 kb	500 µl	1,000
BR332-500	100bp +3K DNA RTU Ladder, 100 bp ~ 3 kb	500 µl	1,000
BR341-500	1 kb DNA RTU Ladder, 250 bp ~ 10 kb	500 µl	1,000
BR343-500	1 kb plus DNA RTU Ladder, 100 bp ~ 10 kb	500 µl	1,000
BR0671	Pre-stained Protein Marker, 9~180 kDa	250 µl x 2	2,400
BR1811	Pre-stained Protein Marker, 9~245 kDa	250 µl x 2	3,000

實驗常用 Buffer

貨號	品名	包裝	價格
BR110-1L	10X PBS Buffer, pH7.4	1L	800
BR130-1L	10X Tris Glycine buffer (10X TG Buffer)	1L	800
BR180-1L	0.5M EDTA	1L	800
BR210-1L	0.5M Tris-HCl, pH6.8	1L	1,000
BR220-1L	1M Tris-HCl, pH6.8	1L	1,000
BR230-1L	1M Tris-HCl, pH7.2	1L	1,000
BR240-1L	1M Tris HCl, pH7.4	1L	1,000
BR250-1L	1M Tris-HCl, pH7.5	1L	1,000
BR260-1L	1M Tris-HCl, pH8.0	1L	1,000
BR270-1L	1.5M Tris-HCl, pH8.8	1L	1,000
BR280-1L	2M Tris-HCl, pH8.0	1L	1,000
BR290-1L	RBC Lysis Buffer	1L	1,500
BR320-200/ 1L	10% SDS Buffer	100ml x2/ 1L	1,000/ 3,000

實驗常用藥品

貨號	品名	包裝	價格
BR420-1KG	Boric acid	1kg	2,800
BR450-1KG	Glycine	1kg	1,800
BR460-1KG	Sodium chloride	1kg	1,000
BR470-1KG	Tris	1kg	2,500



瑞柏生物科技股份有限公司


台北 02-28118200

新竹 03-6660116


台中 04-23152922

高雄 07-5500680

DNA Ladder (核酸標準液)電泳圖

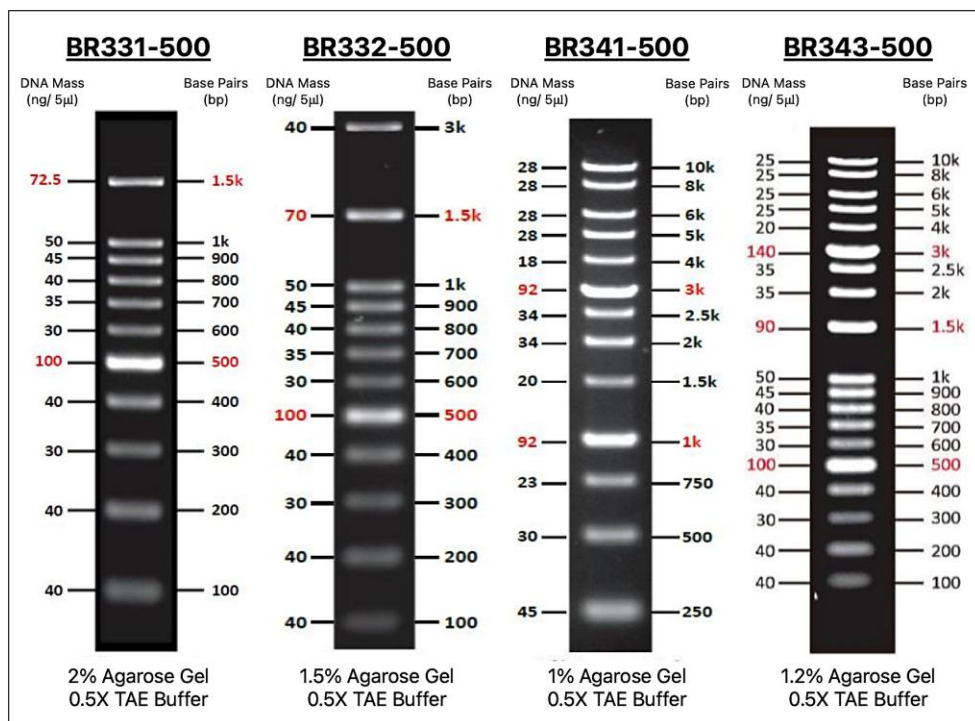

Pre-stained Protein Marker (預染色蛋白質標準品)電泳圖

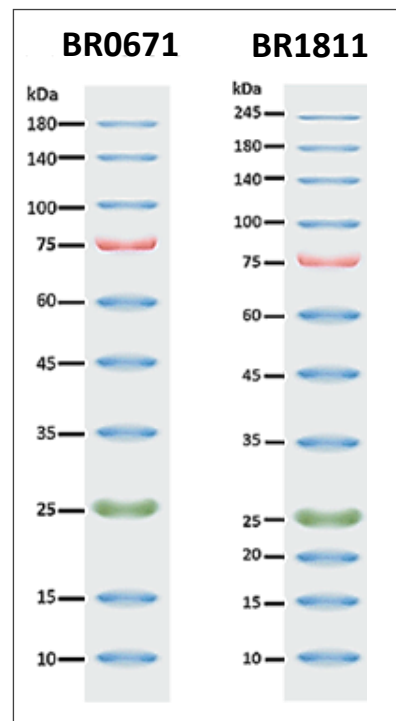


Low EEO agarose (低電滲瓊脂膠)

- **特色**
- 低 EEO (electroendosmosis <0.13)，電泳遷移率高
- 凝膠強度較強 (1% Gel: 1,365 g/cm²)，操作容易、減少膠體破裂
- 解析度較佳，且批次間解析度一致
- 分離效果較好，DNA條帶較平整
- 為 DNase-Free 及 RNase-Free

💡 跑膠時，該選擇 TAE 還是 TBE buffer?

	TAE Buffer	TBE Buffer
成分差異	Acetic acid	Boric acid
導電度	較低	較高
是否易過熱	是	否
可重覆使用	不建議	可
DNA回收	可	不建議
	硼酸會抑制Mg ²⁺ 影響後續實驗	
適用DNA片段	>1kb	<300bp
	介於300bp~1kb之間，則兩者皆可	
二級結構核酸解析度	較好	較差

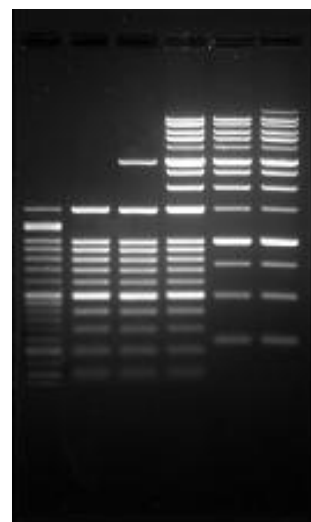

使用 1% agarose
及 1% TAE buffer
跑膠之結果



瑞柏生物科技股份有限公司

台北 02-28118200

台中 04-23152922

新竹 03-6660116

高雄 07-5500680